



AD

AD-E400 169


**TECHNICAL REPORT ARLCD-TR-77031** 

ORGANIC CARBON TESTS ON WASTEWATERS
FROM AMMUNITION PLANTS

PART I. COD AND TOC OF PURE EXPLOSIVES

AD NO.

MILTON ROTH
JOSEPH M. MURPHY, JR.

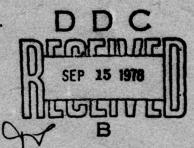

**JUNE 1978** 



US ARMY ARMAMENT RESEARCH AND DEVELOPMENT COMMAND
LARGE CALIBER
WEAPON SYSTEMS LABORATORY
DOVER, NEW JERSEY

APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED.

78 07 26 017

The view, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy or decision, unless so designated by other documentation.

Destroy this report when no longer needed. Do not return it to the originator.

Unclassified
SECURITY CLASSIFICATION OF THIS PAGE (When Date Entered) REPORT DOCUMENTATION PAGE

READ INSTRUCTIONS BEFORE COMPLETING FORM RECIPIENT'S CATALOG NUMBER Technical Report ARLCD-TR-77631 Correlation of Oxygen Demand and Total Organic Carbon Tests on Wastewaters from Final Ammunition Plants, 6. PERFORMING ORG. REPORT NUMBER Part I. COD and TOC of Pure Explosives, Environmental Qual R&D Funds B. CONTRACT OR GRANT NUMBER(+) Milton Roth - ARRADCOM, Project Leader Joseph M. Murphy, Jr - ARRADCOM PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS STB-MTD-LCWSL US Army Armament Research & Development Command Dover, NJ 07801 11. CONTROLLING OFFICE NAME AND ADDRESS JUNE 1078 US ARRADCOM ATTN: DRDAR-TSS Dover NJ 07801
14. MONITORING AGENCY NAME & ADDRESS(II different from Controlling Office) 15. SECURITY CLASS. (of this report) STB-MTD-LCWSL Unclassified US Army Armament Research & Development Command 15a. DECLASSIFICATION/DOWNGRADING Dover, NJ 07801 16. DISTRIBUTION STATEMENT (of this Rep Approved for public release; distribution unlimited. 17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) 18. SUPPLEMENTARY NO 19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Correlation Chemical Oxygen Demand (COD) Explosive Manufacturing Sensitivity Total Organic Carbon (TOC) Precision Statistical Analysis Monitoring Total Organic Nitrogen 28. ABSTRACT (Continue on reverse side if necessary and identify by block number) See reverse side

DD 1 JAN 73 1473 EDITION OF 1 NOV 65 IS OBSOLETE

Unclassified

SECURITY CLASSIFICATION OF THIS PAGE (When Data Bri

410163

# **ABSTRACT**

New Federal regulations for Effluent Limitations, Guidelines, and Standards for Explosives Manufacturing, Point Source Category, has led to the evaluation of the specified oxygen demand tests in order to determine their applicability to munition unique wastewaters. In part one of this study, results of chemical oxygen demand (COD) and total organic carbon (TOC) tests performed on solutions of pure energetic materials are reported and discussed.

The experimental data was analyzed statistically to determine the correlation between the COD and TOC tests when measuring known solutions of energetic materials; the relationship between actual and theoretical values was also determined.

Experimental results show that response to the COD test varies significantly between compounds and is not an accurate indicator of organic loading of munition wastewaters despite its high level of precision. TOC is found to be more accurate and more sensitive.

Unclassified

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

# **ACKNOWLEDGMENTS**

Mr. Nathaniel S. Gelber, Analytical Chemist, U. S. Armament Research and Development Command, Dover, New Jersey provided invaluable assistance by preparing all the standard solutions and performing all the chemical oxygen demand tests. His comments on the details of the testing and his evaluation of the results were most helpful.

This study was funded under the Environmental Quality R&D program administered by the DARCOM Lead Laboratory, Aberdeen Proving Ground, MD.

ACCESSION f	or
NTIS	White Section
DDC	Buff Section
UNANNOUNCE	ED 🗆
JUSTIFICATIO	
DISTRIBUTION	N/AVAILABILITY CODES
Dist. AVA	IL. and/or SPECIAL

# TABLE OF CONTENTS

Page	NO
Introduction	1
Experimental	2
Experimental Design Procedure	2 2
Discussion of Results	3
Comparison of COD and TOC for Nitrotoluenes and Toluene at High Concentrations	3
Comparison of COD and TOC for Nitrotoluenes and Toluene at Low Concentrations	3
Comparison of COD and TOC for Nitrobenzenes and Benzene at High Concentration	4
Comparison of COD and TOC for Nitrobenzenes and Benzene at Low Concentration	4
Comparison of COD and TOC Tests for Nitramines, Guanidine Nitrate, and Urea at High Concentration	4
Comparison of COD and TOC Tests for Nitramines, Guanidine Nitrate, and Urea at Low Concentration	5
Comparison of COD and TOC for Nitrate Esters, Alcohols, and Acids at High Concentration	5
Comparison of COD and TOC for Nitrate Esters, Alcohols, and Acids at Low Concentration	5
Experimental Data	5
Conclusions	6
Recommendations	7
References	7
Distribution List	35

# Tables

1	Point sources from manufacturing and loading of explosive	es 9
2	EPA effluent limitations, guidelines and standards	10
3	COD and TOC test results for toluene and nitrotoluenes	11
4	COD and TOC test results for benzene and nitrobenzenes	16
5	COD and TOC test results for urea, nitramines, and guanidine nitrate	18
6	COD and TOC test results for nitrate esters	21
7	COD and TOC test results for alcohols and acids	22
8	Statistical parameters for aromatic compounds	24
9	Statistical parameters for nonaromatic compounds	25
10	Comparison of completeness of recovery by COD and TOC tests	26
'igui	res	
1	Comparison of COD and TOC for nitrotoluenes and toluene at high concentration	27
2	Comparison of COD and TOC for nitrotoluenes and toluene at low concentration	28
3	Comparison of COD and TOC for nitrobenzenes and benzene at high concentration	29
4	Comparison of COD and TOC for nitrobenzenes and benzene at low concentration	30
5	Comparison of COD and TOC tests for nitramines, guanidine nitrate, and urea at high concentration	31
6	Comparison of COD and TOC tests for nitramines,	32

7	Comparison of COD and TOC for nitrate esters, alcohols and acids at high concentration	33
8	Comparison of COD and TOC for nitrate esters, alcohols and acids at low concentration	34

### INTRODUCTION

Effluent limitations and guidelines for existing sources, to be achieved by the application of the best practicable control technology currently available (BPCTCA), and the best available technology economically achievable (BATEA), were promulgated by the Environmental Protection Agency (EPA) in the Federal Register of 9 March 1976 for Explosive Manufacturing Point Source category. These regulations have led to an evaluation of the test methods for determining organic loading and the resulting oxygen demand of military-unique wastewaters. The interim final rule for effluent guidelines, for explosive manufacturing point source guidelines sets standards in terms of the following parameters: biochemical oxygen demand (BOD<sub>5</sub>), chemical oxygen demand (COD), and total suspended solids (TSS).

The energetic materials found in wastewaters of different types of explosive manufacturing are listed in table 1. In the first part of this study, known concentrations were tested for COD and total organic carbon (TOC). The values obtained were statistically analyzed to determine which method is more appropriate for monitoring waste streams containing energetic materials and also to determine if any relationship exists between the methods.

The discharge standards set for explosive manufacturing point source category are listed in table 2 (ref. 1). These standards, in terms of kilograms (lb) BOD or kilograms (lb) COD allowed per 1000 kg (2200 lb) of product produced are established for BPCTCA and for BATEA. When the proposed discharge standards are projected over maximum production rates for nitroglycerin and trinitrotoluene and the water usage associated with their production, (ref. 2) the resulting discharge concentrations, in milligrams per liter  $(mg/\ell)$ , of COD in the waste streams are extremely low (table 2). These standards present a problem of how to measure these parameters at such low levels. The COD analysis has limited accuracy below 10 mg/l of COD (ref.3). In order to evaluate the EPA environmental water quality standards for munition wastewaters, the applicability of the conventional oxygen demand parameters to low concentration levels of COD of military-unique compounds was studied. In this study, consideration was given to other tests capable of detecting these munition-unique pollutants. Thus, the TOC test is compared to the COD for the same compounds at the same concentration.

An evaluation of this type will not only show the precision and accuracy of each method, but also will indicate which is the most appropriate for monitoring explosive manufacturing point sources.

## **EXPERIMENTAL**

# Experimental Design

Since one objective of this study was to compare different test methods to evaluate their suitability for estimating compliance with the EPA standards, a factorial design was used with the intention of applying a variety of statistical analyses. In this design, the following factors were considered:

- 1. Method (COD and TOC)
- 2. Concentrations (high and low levels)

Four replications were conducted on each combination of factors. The replications were randomized to eliminate the effect of time on the results. Analyst effect was eliminated by having one analyst perform all the COD tests and a second analyst performed all the TOC tests. The test solutions were prepared as described in the procedure paragraph below and aliquots were taken from the same solution for all tests. In general, the COD tests were completed before beginning the TOC tests.

Prior to conducting the analysis of variance of the data, the variances were checked for homogeneity by the Bartlett test (ref. 5). The result showed that the variances were not homogeneous; therefore, the analysis of variance, which is based on the assumption of homogeneity, was not performed. The statistical tests that were performed are discussed later in this report.

### Procedure

Standard solutions that are typically found in the wastewaters of army ammunition plants were prepared. Isomers of the compounds were also included in the experiment in order to examine the effect of the structure of a compound on the test response. Each standard solution was prepared by dissolving weighed amounts (to the nearest 0.1 mg) of the compounds in distilled water. The mixture was stirred with a magnetic stirrer until solution was complete. Several compounds required more than two days of constant stirring to effect solution.

The solutions were cooled to ambient temperature before diluting to the mark with water. Aliquots were transferred from the high concentration level solutions and diluted appropriately to prepare the low concentration level.

The chemical oxygen demand analyses were performed according to the standard method (ref. 6). The total organic carbon analyses were also performed according to the standard method (ref. 7) using a Beckman 915-A Total Organic Carbon Analyzer, a Hamilton spring-loaded syringe, and a sample size of 20 micro-liters.

Theoretical COD values were used as a basis for calculating the percentage of recovery for the different compounds. The theoretical COD values were calculated on the assumption of complete oxidation of all elements in the energetic materials as follows:

$$\begin{array}{cccc} 2H & \longrightarrow & H_2O \\ C & \longrightarrow & CO_2 \\ 2N & \longrightarrow & N_2O_5 \end{array}$$

Theoretical TOC values were used as the basis for calculating the percent of recovery. The percentage values were calculated on the assumption of complete combustion of the carbon in the compound to carbon dioxide.

### **DISCUSSION OF RESULTS**

Comparison of COD and TOC for Nitrotoluenes and Toluene at High Concentrations

As shown in figure 1, at the high concentration level the TOC concentration shows a recovery of approximately one hundred percent. Toluene yielded low results, probably because of evaporation during the sample preparation (ref. 8).

The COD test at the high concentration level displayed results that varied with each compound. Low results were displayed by m-nitrotoluene; 2, 3 dinitrotoluene; 3, 4 dinitrotoluene; and toluene. Note that the nitrotoluenes yielding low results all had a nitro group in the meta position. All the other results were  $\pm$  5% from 100% recovery.

Comparison of COD and TOC for Nitrotoluenes and Toluene at Low Concentrations

At low concentrations (fig. 2), the mononitrotoluenes all yielded recoveries that range from 95% - 105% in both TOC and COD. These results indicate that the mononitrotoluenes can be analyzed with good recovery by either test method. However, the more highly nitrated toluenes were significantly lower by the COD test than by the TOC.

An analysis of a solution containing 6.41 ppm TNT (equivalent to 2.22 mg/L TOC) is likely to produce erratic results since the detection limit, under ordinary procedures, is limited to 1 mg/L (ppm) of TOC (ref. 9).

# Comparison of COD and TOC for Nitrobenzenes and Benzene at High Concentration

At the high level (fig. 3), TOC and COD deviate significantly from full recovery for most of the compounds tested. The COD results of this group are consistently low, with the dinitrobenzenes registering very minimal COD values. Nitrobenzene and benzene also yielded low results. This could be the result of the refractoriness of aromatics which typically are not oxidized appreciably in COD analysis.

# Comparison of COD and TOC for Nitrobenzenes and Benzene at Low Concentration

At the low level concentration (fig. 4), good recoveries for TOC are exhibited by the nitrobenzenes. Benzene again yielded a low result, probably due to losses from volatilization of the compound.

The COD results follow the same pattern as at the high level, except that the nitrobenzene recovery is now complete. This indicates that either the concentration of dichromate or the digestion time at the high level is inadequate.

# Comparison of COD and TOC Tests for Nitramines, Guanidine Nitrate, and Urea at High Concentration

All the compounds tested by the TOC method yielded low results except for tetryl and guanidine nitrate (fig. 5). The COD test also displayed consistently low results for this group of compounds. Nitroguanidine, guanidine nitrate, and urea did not respond at all to the COD. Although it has been recognized that certain aromatics give low results in the COD test, this appears to be the first time it has been found that the urea structure resists dichromate oxidation.

Tetryl and RDX displayed low COD results which indicates that nitramines also respond incompletely to the COD test. Tetryl is higher than the rest, probably because it contains three nitro groups in addition to a nitramine.

# Comparison of COD and TOC Tests for Nitramines, Guanidine Nitrate, and Urea at Low Concentration

TOC recoveries were satisfactory for HMX and RDX, but high results were obtained on nitroguanidine and guanidine nitrate, and low results were obtained on tetryl and urea (fig. 6). The COD results followed the same pattern found at the high concentration level, indicating that the effect of structure is independent of concentration for this group of compounds.

# Comparison of COD and TOC for Nitrate Esters, Alcohols, and Acids at High Concentration

Both TOC and COD (fig. 7) appear to be very appropriate for monitoring these compounds. Glycerin, pentaerythritol, and acetic acid are more typical of the type of compounds found in industrial rather than explosives wastewaters and are at least part of the basis for assuming that COD is a representative measure of organic loading in a stream. However, the data from explosives and urea type compounds indicates that the COD test has still-to-be-defined limitations.

# Comparison of COD and TOC for Nitrate Esters, Alcohols, and Acids at Low Concentration

TOC results (fig. 8) are erratic, probably because the measurements were made at the lower level of sensitivity of the TOC analyzer. COD results, with the exception of acetic acid, yielded good recoveries with these more conventional materials. Acetic acid yielded a low recovery, probably because of losses due to volatilization during the dichromate oxidation.

## **Experimental Data**

The experimental data is reported in tables 3 through 7. The tables include the standard deviations and the percent recoveries of COD and TOC test results. In order to complete the study it was necessary to determine the correlation between the test methods to see if the methods could be interchanged. Correlation coefficients were calculated using the recovery values obtained from the COD and TOC test methods of the combined high and low levels of concentration for individual compounds. Significant correlations at the 95% confidence level were obtained by several compounds in each table, but no overall correlation was found. Results from the statistical analyses are in tables 8 and 9.

Generally, in both the aromatic and nonaromatic groups, variances were always higher on results from total organic carbon analyses as opposed to chemical oxygen demand analyses. The larger variances might be an indication of the greater sensitivity of the TOC test method and the need for replication in order to obtain a reliable estimate of the true value.

As a firal test of the overall completeness of recovery of the COD and TOC test methods, t-tests were performed to compare the data from each test with the theoretical recovery value of one hundred percent. The statistical evaluation (see table 10) showed that there was a significant difference from 100% recovery at the 90% level of confidence with 23 degrees of freedom for the COD analyses at both high and low concentration levels. With the TOC analyses, however, there is no significant difference from 100% recovery at either concentration level. This indicates that the TOC test is more accurate and appears to be a more appropriate measure of the organic loading of a wastewater.

### CONCLUSIONS

- 1. The results obtained by the COD test do not give complete recovery with many compounds. Although this fact is recognized for aromatic compounds and is confirmed in this study, the study also revealed that compounds with urea structure, nitramine groups and certain nitro groups yield low results. However, the precision of the COD method surpasses that of the TOC.
- 2. The results obtained by the TOC test do give complete recovery as indicated by the t-test. The principle of the TOC test makes it independent of the structure of a compound. However, the precision of the TOC test, particularly at low concentrations, would require numerous replications to obtain a reliable estimate of the true mean of a sample.
- 3. From the point of view of enforcement of standards, the COD test appears inappropriate for the explosives industry. The TOC test would be more meaningful from the point of view of accuracy, but further study of this test is needed to improve its precision.
- 4. Although it has been shown that TOC and COD tests are correlated for many industrial wastes (ref. 10), their correlation is not a general rule because of the widely different principles upon which the tests are based. The data presented in this study indicates that correlations occur by chance rather than by a unifying principle. Thus, the existing standard based on the COD test for the explosive industry, cannot be converted to an equivalent standard based on the TOC test.

### RECOMMENDATIONS

- 1. Further study should be made of the COD test to determine the effect of molecular structure on oxidizability. The efficiency of other oxidants, as compared to dichromate, should also be studied.
- 2. Consideration should be given to replacing the COD test with TOC in the EPA Limitations, Guidelines and Standards for the Explosives Industry.
- 3. Further study of the general indicators of organic loading in wastewaters and receiving waters is needed. For the explosives industry, a combination of TOC and total organic nitrogen should be considered to obtain a more pertinent measure of water quality.
- 4. Actual samples of wastewaters from explosive plants should be tested to confirm results obtained from the theoretical studies.

### REFERENCES

- Federal Register, Vol 41, No. 47 Tuesday, March 9, 1976, pp. 10180 - 10188
- D. H. Maybury, J. L. Evans, and A. Carotti, "Propellant Plant Pollution Abatement Improvement of Water Utilization at Radford Army Ammunition Plant," Picatinny Arsenal Report 4562
- "Standard Methods for the Examination of Water and Wastewater,"
   14th Edition, American Public Health Association, Washington, DC,
   1976, p. 551
- J. P. Glennon, CPT, MSC, L. Reuter, LTC, MSC, "Environmental Quality Standards for Munitions - Unique Pollutants," U. S. Army Bioengineering Research and Development Laboratory, Fort Detrick, Frederick, Maryland. Presented at Seventeenth DOD Explosives Safety Seminar, Denver, Colorado, September 14, 1976
- C. A. Bennett and N. L. Franklin, "Statistical Analysis in Chemistry and the Chemical Industry," Wiley and Sons, Inc., NY, 1954, p. 197
- "Standard Methods for the Examination of Water and Wastewater,"
   14th Edition, American Public Health Association, Washington, DC,
   1976, pp. 550 555

- 7. Ibid, pp. 532-534
- C. E. Van Hall, D. Barth, and V. A. Stenger, Analytical Chemistry, Vol 37, No. 6, May 1965, pp. 769-773
- "Standard Methods for the Examination of Water and Wastewater,"
   14th Edition, American Public Health Association, Washington, DC,
   1976, p. 533
- 10. Environmental Protection Agency, "Handbook for Monitoring Industrial Wastewater," August, 1973, pp. 5-10

Table 1. Point sources from manufacturing and loading of explosives

# TNT Production

3,4-INT, 3,5-DNT, mononitrotoluene, 1,3-dinitrobenzene Condensate Water: 2,3-dinitrotoluene (DNT), 2,4-DNT, 2,5-DNT, 2,6-DNT,

# TWT Load, Assembly and Pack

Pink Water: 2,4,6-trinitrotoluene (TNT), 2,6-DNT, 4,6-DNT

# RDX/HMX

Wastewaters: Cyclotrimethylene trinitramine (RDX), cyclotetramethylene tetranitramine (HMX)

# Mitrate Ester Wastewaters

Trinitroglycerol (NG), Triethyleneglycoldinitrate (TEGDN)

Table 2. EPA effluent limitations, guidelines and standards\*

	LIMITATIO	NS, kg/kkg C	LIMITATIONS, kg/kkg OR LBS/1000 LBS PRODUCT	RODUCT
EFFLUENT CHARACTERISTIC	FOR ANY 1 DAY, MAX BPCTCA BATEA	Y, MAX BATEA	AV OF 30 CONSEC DAYS, MAX BPCTCA BATEA	C DAYS, MAX BATEA
COD	1	0.85	•	0.55 b
<b>EXISTING SOURCES</b>	17.7		2.59	
NEW SOURCES	3.6		2.3	
BOD5	1	0.11	•	0.067
<b>EXISTING SOURCES</b>	0.72		0.24	
NEW SOURCES	0.35		0.23	
COD @ MA)	COD @ MAX PROD. RATES:	NG, mg/l TNT, mg/l	a b 26 17 5.4 3.5	

\* 9 MAR 1976

Table 3. COD and TOC test results for toluene and nitrotoluenes

Commodiand	Concentration	COD.	1/200	Concentration	TOC, II	1/8	Recove	ry, \$
Component	mg/1	Added Found	Found	mg/l	Added Found	Pound	8	1300
Toluene	169.5	529.9	408.0 385.0 391.0 362.0	169.5	154.5	54.0 33.8 17.2 33.8	76.9 72.6 73.7 68.3	31.8 19.9 10.1
		<b>×</b> @	387.0		<b>(</b> × ∞	34.7	3.60	20.4 8.88
Toluene	17.0	53.0	31.1 33.0 28.7 29.4	17.0	15.5	7.9 8.0 9.7	58.6 82.6 55.5 5.5	51.0 51.7 51.0
		<b> </b>   0	30.6		p< ∞	7.7	3.6	3.03
S. S	204.2	428.8	418.4 418.4 418.4 419.4	200.6	123.0	120.0 122.0 123.6 118.0	94.5 94.5 94.8	98.5 186.0 95.9
		l× ø	418.7		ا≿ ∞	120.9	97.6	98.2
TH-0	20.4	42.9	42.1 39.7 40.1	20.1	12.3	11.5 13.5 14.8	88.88.98 6.54.63	93.4 109.7 105.6 120.3
		<b>l</b> ≻ ∞	1.2		j× ¤	13.2	95.3	107.3

Table 3. (Continued)

Compound Concentration	1/8m	п <b>-ит</b> 192.9		п-ит 19.3		p-NT 31.6		2-NT 6.10	
COD, mg/1		405.1	l× æ	40.5	l× ∞	ħ*99	<b>×</b> &	12.8 10.5 11.4	ı×
mg/1	round	336.9 342.9 347.0 341.0	342.0	37.4 38.5 38.4	38.5	62.7 63.1 64.7	63.5	12.6	7.11
Concentration	mg/1	193.8		19.4		31.6		6.34	
TOC, mg/1	Added	118.8	<b>&gt;</b> 0	11.9	l≍ ø	19.3	<b> </b>   0	3.89	ı×
mg/1	Lound	134.0 131.6 136.0	134.4	11.0	3.0	19.0 19.0 19.5	19.3	3.6 3.0 3.0	3.6
Recov				8484 				888.48 4.0.2.0	
Recovery, \$	3	112.7	113.1	8888 26.20	93.2	98.7 98.7 101.3 101.3	100.0	165.3 98.5 95.1	8.5

Table 3. (Continued)

Compound	Concentration	COD.	ma/1	Concentration	TOC. II	1/2	Recov	erv. \$
	mg/1	Added	ded Found	mg/1	Added Found	Found	8	TOC
2,3 DNT	50.5	79.8	62.8 63.2 62,7 57.8	50.5	23.4	25.98 26.98 26.08	78.7 78.6 72.4	103.8 105.7 105.7
		j∺ ∞	61.6		<b> </b> × ∞	24.8	3.23	105.8
2,3 DNT	6.38	10.1	5.5 4.0 7.0	6.38	2.95	1.9	66.3 53.5 61.4 69.3	67.8 57.6 64.4 57.6
		× @	6.3		<b>K</b> #	1.83	62.6	62.0 5.1
2,4 INT	6.66	157.9	153.0 154.0 156.0 157.0	6.66	46.3	49.0 49.0 48.8 48.0	96.8 97.5 99.4	105.8 105.8 105.3 103.6
		l∺ æ	155.0		<b>P</b>	148.7	98.1	105.2
2,4 DNT	6.7	10.6	10.2 10.1 8.9 8.7	6.7	3.10	8.6.6.6	883.58 8.6.6.6	88.00 87.0 9.00
		<b> </b>	9.0		<b> </b> ≻ =	2.8	4.68	8.4

Table 3. (Continued)

Como	not tent mooned	900		Contract to		<b>'</b>	December	,
o de la composición della comp	mg/1	Added Found	Found	mg/1	Added Fo	Found	COD	TOC
2,5 DNT	50.5	79.8	75.5 77.6 75.2 4.6	50.5	23.4	25.0 24.5 24.5 2.45	8298 3.00.00	106.8 104.7 104.7
		j× ∞	76.2		<b>×</b> ∞	24.6	95.5	105.3
2,5 DNT	6.45	10.2	8.9 9.4.9 4.6.4.9	6.45	2.98	3.0	88.5.5.5 84.3.5.5	83.8 100.7 134.2 100.7
		l∺ ø	9.1		<b> ≻</b> ∞	3.13	9.0	21.1
2,6 DNT	50.8	80.3	77.8 79.2 77.9	50.8	23.5	24.7 23.8 24.7	98.98.9	105.1 101.2 109.8 105.3
		<b> </b>   0	78.6		<b> </b>	24.8 0.82	97.9	105.4 3.52
2,6 DNT	84.9	10.2	6.9 9.9 9.9	6.48	3.0	0.4.0.	93.0 97.0 97.0	100.0 123.3 106.7 100.0
		K @	8.5		H @	3.23	83.0	107.5

Table 3. (Continued)

Compound	Concentration mg/1	Added Found	mg/1 Found	Concentration mg/l	Added Found	Found	Reco	Recovery, \$
3,4 DNT	53.8	85.1	56.6 55.9 54.9 53.2	53.9	25.0	%%%% %.0.%%	88.68 6.68 6.59 8.50 8.50 8.50 8.50 8.50 8.50 8.50 8.50	107.4 104.2 102.2 104.2
		H @	55.2		p¥ w	26.1	64.8 1.72	104.5 2.15
3,4 DNT	6.23	8.6	6.54 5.1.69	6.23	2.88	3.00.0	60.22 60.03 60.03	72.9 4.69.4 121.5
		l× ∞	5.3 0.77		<b> </b>	4.0	54.1 7.85	83.3
2,4,6 INT	1.00.4	127.3	124.1 123.1 124.1 120.8	100.4	37.2	37.0 37.0 37.0	2824 ~	8848 2.5.1.2
		p× œ	123.0		p× ∞	35.9	96.6	3.57
2,4,6 THT	6.41	8.13	6.6 6.3 7.7 7.7	6.41	2.28		85.68 5.69	139.6 202.0 90.1 180.1
		<b>×</b> ∞	1.1		<b>p</b> < ∞	3.4	86.7	153.0

Table 4. COD and TOC test results for benzene and nitrobenzenes

Recovery, \$	67.9 63.7 59.0 61 59.0 61		33.0 37.6 33.0 36.1 36.8 23.5 15.4		8.39 111.8 7.50 90.1 5.69 90.1 6.20 90.1			
			44.00 v.v.o.i				0.6.0.0	
TOC, mg/l Added Found	136.2	(H @	13.6	[ <b>1</b> € 6	4.40	l⊷ ø	<b>11.</b> 9	<b>№</b> •
Concentration mg/1	148.0		14.8		150.2		15.02	
COD, mg/l Added Found	307.0 288.0 267.0 270.0	283.0	17.6 17.0 16.3 10.6	15.4	18.0 16.1 12.2 13.3	14.9	0100	1.2
Added	452.0	K @	45.2	<b>≻ </b>	214.5	p< ∞	21.45	p× e
Concentration mg/l	147.1		7.41		150.2		15.02	
Compound				Nitrobenzenes	m-DN Benzene		в-IN Benzene	

Table 4. (Continued)

Compound	Concentration mg/1	COD, mg/l Added Found	mg/1 Found	Concentration mg/1	TOC, mg/1	Found	Mecon COD	Recovery, \$
p-IN Benzene	51.3	13.2	6.6 7.5 6.3	51.3	22.0	45.5 21.6 19.0 24.5	9.0 10.2 8.6	207.0 98.3 86.4 111.5
		[¥ ∞	7.1		<b>×</b> ∞	27.7	9.61	125.8
p-DN Benzene	15.4	22.0	20.00	15.4	9.9	6.8 6.0 6.9	5.4 9.09 12.3 10.9	103.0 90.9 110.6 104.5
		l <b>⋉</b> æ	2.1 0.7		<b>×</b> &	6.75	9.42	102.3
Nitrobenzene	156.7	305.4	198.0 180.0 216.0 201.0	156.7	7.16	83.0 80.0 99.0	64.8 78.9 65.8	90.5 87.2 87.2 107.9
		j⊯ ∞	199.0		l≍ œ	85.5 9.11	65.1 4.84	93.2
Nitrobenzene	15.7	30.5	29.9 31.8 30.5	15."	9.17	0.00 0.00 0.4 1.00	97.9 104.1 97.5 99.9	%%% 6.7%
		l∺ œ	30.5		<b> </b>    0	8.68	3.05	9.4° 2.99

Table 5. COD and TOC test results for urea, nitramines, and guanidine nitrate

Recovery, \$	77.1 79.5 79.5	3.66	52.8 48.0 84.1	61.2	78.0 84.1 8.8	81.0 4.0	72.5 57.9 217.3 144.9	123.0 73.4
Reco	0 -0.36 0.53 0.18	0.09	0.86 0.86 -1.76 -3.58	0.90 2.16	0.13 0.17 0.21 -0.13	0.0	-2.16 -3.6 -1.08 7.38	3.55
mg/1 Found	33.00 33.00 33.00 33.00	31.3	9998 905.	2.55 0.67	45.0 44.5 48.5 49.0	46.8	0.01 2.01 0.01	0.85
TOC, mg/l	41.5	<b>×</b> ∞	4.16	<b> </b> •	57.7	<b>×</b> ∞	69.0	<b>l</b> κ α
Concentration mg/l	207.6		20.8		502.2		6.03	
COD, mg/l	0.00 -1.58 2.35 0.79	0.39	-0.38 -0.38 -0.78 -1.58	-0.78 0.57	1.17	0.88	-0.24 -0.40 -0.12 -0.82	0.02
COD,	7.044	× &	14.1	<b> </b>	927.3	l× ∞	11.11	Þ×∞
Concentration mg/1	206.8		20.7		502.6		. 20.9	
Compound	Urea		Ures		NG.		NG.	

Table 5. (Continued)

ery, \$	71.3 69.3	70.8	99.0 148.5 99.0	24.8	98.2 93.5 106.5 93.5	96.5 47.9	7.4.7.4 7.4.7.4	75.0 22.5
Recovery, \$	8888 8.5.6.6 8.6.6.6	1.60	88.88 83.68 83.68	28.0	78.7 81.9 82.4 81.1	81.1	57.1 64.1 87.1	68.9
Found	33.6	3.6	1.0	1.13	13.8 14.0 16.0 14.0	14.5	0.1.1.	1.35
TOC, mg/1 Added Found	5.05	p× ∞	1.01	<b> </b>	15.0	p× w	1.80	l∺ ø
Concentration mg/l	31.0		6.20		51.1		6.14	
Found	13.2 12.8 11.8 13.0	12.7	2.6 1.9 1.9	1.9	47.1 49.0 49.3 48.5	48.5 1.0	444 6.4 6.5	4.95
COD, mg/l	39.4	<b>×</b> ∞	3.04	× &	59.8	<b> </b>   0	7.18	<b>l</b> ≍ ∞
Concentration mg/1	30.4		6.20		51.1		4۲.9	
Compound	RDX		RDX		Tetryl		Tetryl	

Table 5. (Continued)

		200	,			4		,	
Compound	Concentration mg/l	Added Added	Added Found	concentration mg/l	Added R	Round	COD	OD TOC	
них	5.8	7.58	7.0 6.8 7.5	5.8	₽.0	1.0	92.3 81.8 89.7 98.9	53.2 106.3 106.3 138.3	
		× @	6.88		l× ∞	0.95	7.07	101.2 35.3	
Guanidine Nitrate	300.9	473.0	0.08	301.0	29.6	28.8 30.5 30.0 29.8	0.17 -0.17 -0.25	97.3 103.0 10013	
		<b> </b>	-0.3		× &	29.8	~0.06 0.19	100.6	
Guanidine Nitrate	30.12	4.74	4.4.0 000	30.1	3.47	5.0.0 5.0.0 5.0.0	8.0°0°	149.8 144.0 146.9	
		K @	0.1		<b>×</b> ∞	5.08	0.80	146.18 2.78	

Table 6. COD and TOC test results for nitrate esters

Recovery, \$	105.6 100.6 103.1 89.1	99.7	106.0 127.3 127.3 159.0	130.0	128.3.3 2.5.5.7.	103.3	4.488.55 4.4.98.55	103.3
Reco	98.5 98.5 99.5 99.5	98.7	140.4 75.3 92.5 106.2	103.6 27.6	100.5 99.2 98.8 99.3	99.4	8428 4.2.0.0.	96.6
g/1 Found	\$884 0.000	79.3	1.00	1.23	135.0 135.0 141.0 180.0	148.0	13.5 13.5 14.1 18.0	14.8 2.17
TOC, mg/l Added Found	79.5	j⊭ œ	46.0	<b>×</b> ∞	143.2	K @	14.3	K @
Concentration mg/l	500.0		5.93		479.0		6.74	
ng/1 Found	242.8 242.9 242.9 241.6	243.3 0.9	8.1 3.7 3.1	0.80	481.3 475.2 473.0 475.4	176.2 12.6	48.8 45.3 45.5	46.4 1.3
COD, mg/l	246.6	× &	2.92	× 0	478.9	j× ∞	47.93	<b> </b>
Concentration mg/l	500.0		5.93		479.2		47.96	
Compound	NG		NG		TECTM		TECON	

Table 7. COD and TOC test results for alcohols and acids

Recovery, \$ 100	92.7 92.5 92.5	94.2	8888 6.6.6.6	94.6	99.0 100.7 96.8 100.6	99.3	81.2 92.3 114.4 136.5	24.5
COD	99.0 99.0 9.86	99.1	96.9 100.6 97.9 95.9	97.9	98.0 98.0 97.8 98.2	98.0	85.5 117.7 86.6 98.2	97.0
mg/1 Found	145.3 150.0 145.0 150.0	147.5 2.80	14.5 14.5 15.0 15.0	14.8	132.0 134.2 129.0 134.0	132.3	9.98.88 9.45.45	2.88
TOC,	156.7	× ∞	15.6	l≍ ∞	133.21	<b> </b>	<b>4.</b> 8	× &
Concentration mg/l	398.0		39.8		302.0		6.14	
mg/l Found	477.1 476.4 478.5 476.1	477.0	46.7 48.5 47.2 46.4	47.2	417.0 417.0 416.0 418.0	417.0 0.82	7.4 10.2 7.5 8.5	8.4 1.3
Added	481.4	<b>×</b> ∞	7°84	<b> </b>	425.5	l× w	8.66	l <b>≻</b> ∞
Concentration mg/l	395.9		39.6		301.7		6.14	
Compound	Glycerine		Glycerine		Pentaerythri tol		Pentaerythri tol	

Table 7. (Continued)

Compound	Concentration mg/1	Added	mg/1 Found	Concentration mg/1	TOC,	mg/1 Found		ery, \$
Acetic Acid	313.2	333.9	328.0 326.0 325.0 327.0	313.2	125.3	120.0 123.6 120.0 120.0		95.8 98.6 95.7 7.7
		l <b>≍</b> ∞	327.0		<b>×</b> ∞	120.9		96.5 1.44
Acetic Acid	69.9	7.13	744 E	69.9	2.68	1.00.1		71.5 78.4 66.3 59.7 56.3 37.3 54.6 37.3
		j× ∞	4.4		<b>₩</b> @	1.43	-	53.2

Table 8. Statistical parameters for aromatic compounds

Correlation Coefficient	\$6.00.00.00.00.00.00.00.00.00.00.00.00.00
Variance COD TOC	77.1 282.2 4.75 78.1 33.9 116.6 28.3 75.2 85.7 75.2 85.7 75.2 193.2 72.6 21.9 191.5 193.2 60.5 60.5 409.7 61.6 1956.3 282.2 234.1 5.86 97.4 4.37 1490.0 359.9 46.6
Compound	Toluene o-Nitrotoluene m-Nitrotoluene 2,3 Mnitrotoluene 2,4 Mnitrotoluene 2,5 Mnitrotoluene 2,6 Mnitrotoluene 2,6 Mnitrotoluene 2,4,6 Trinitrotoluene Benzene m-Unitrobenzene p- Unitrobenzene p- Unitrobenzene Nitrobenzene Nitrobenzene

\* Values significant at 95% confidence level

Table 9. Statistical parameters for nonaromatic compounds

Correlation Coefficient	-0.48 0.35 0.30 -0.72 -0.145 0.055 *
TOC	171.6 2819.6 599.8 734.4 1246.1 196.0 2.99 272.3
Variance COD	2.34 10.4 0.31 17.2 50.0 334.9 4.67 8.25 96.0
Compound	Urea Nitroguanidine Guanidine Nitrate RDX HMX Nitroglycerin Triethylene Glycol Dinitrate Glycerin Pentaerythritol Acetic Acid

\* Values significant at 95% confidence level

Table 10. Comparison of completeness of recovery by COD and TOC tests

t - TEST VALUES

TOCLO	48.0
TOC <sub>H1</sub>	1.55
$^{\text{OI}}_{\text{IO}}$	5.76*
COD <sub>H1</sub>	¥99° 1

\* VALUES ARE SIGNIFICANTLY DIFFERENT FROM 100% RECOVERY AT 90% LEVEL OF CONFIDENCE FOR 23 DEGREES OF FREEDOM.

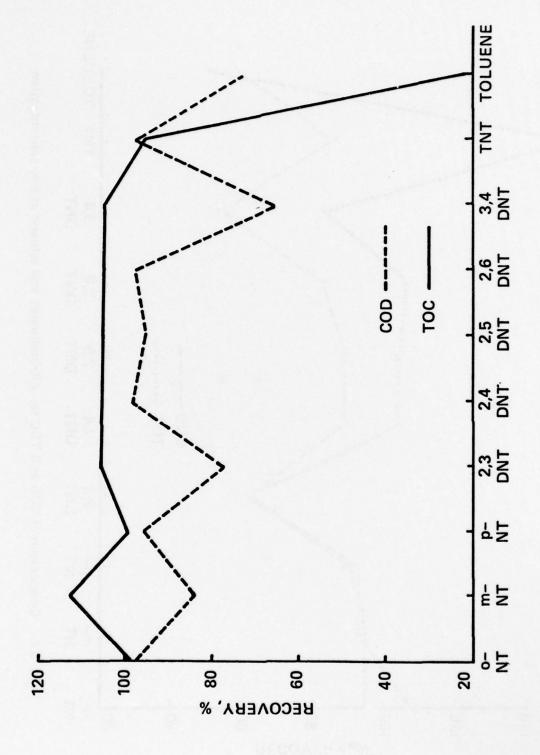


Figure 1. Comparison of COD and TOC for nitrotoluenes and toluene at high concentration

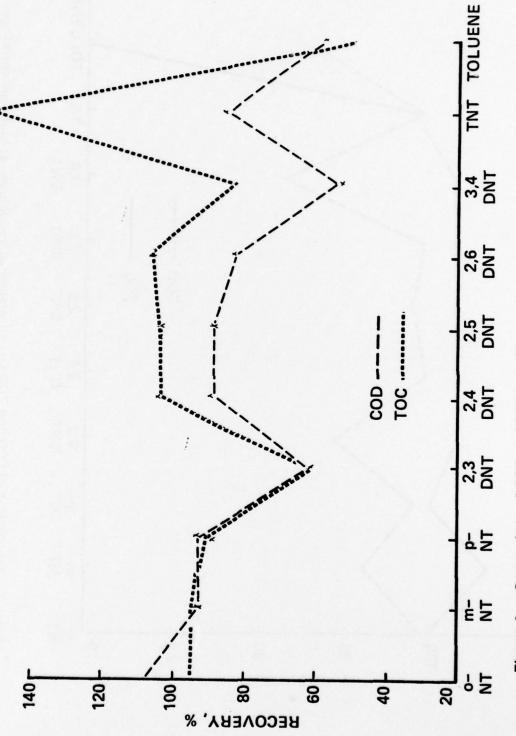


Figure 2. Comparison of COD and TOC for nitrotoluenes and toluene at low concentration

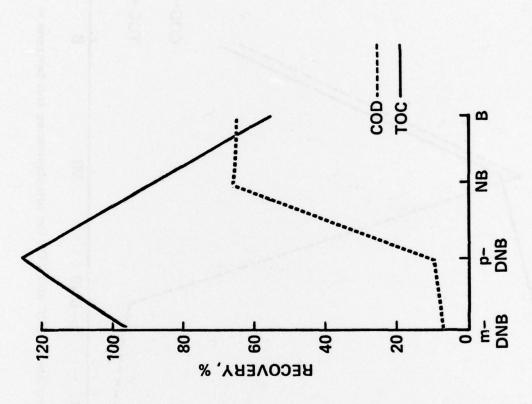


Figure 3. Comparison of COD and TOC for nitrobenzenes and benzene at high concentration

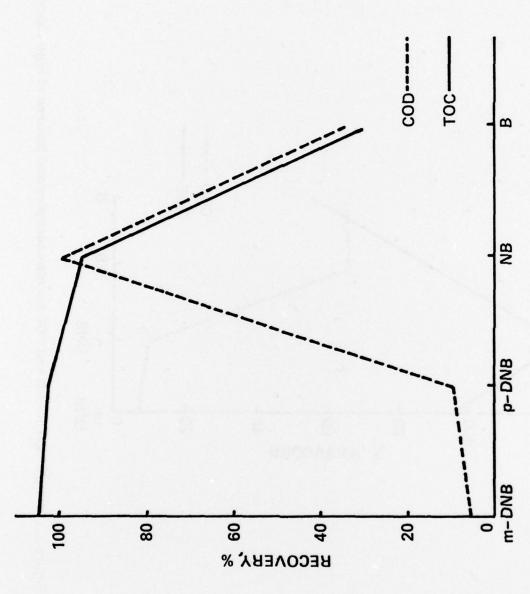


Figure 4. Comparison of COD and TOC for nitrobenzenes and benzene at low concentration

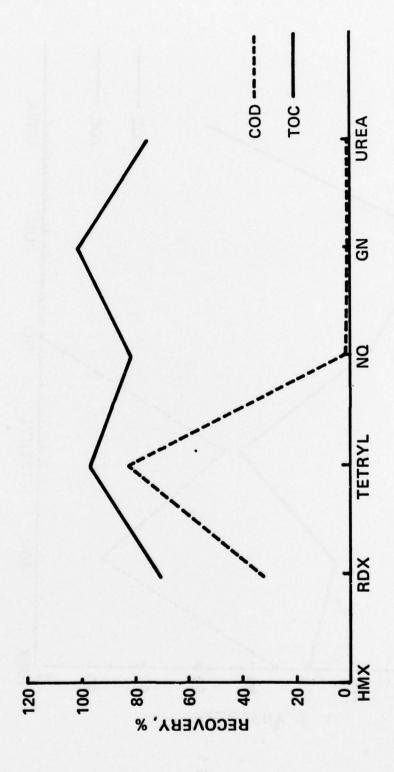


Figure 5. Comparison of COD and TOC tests for nitramines, guanidine nitrate, and urea at high concentration



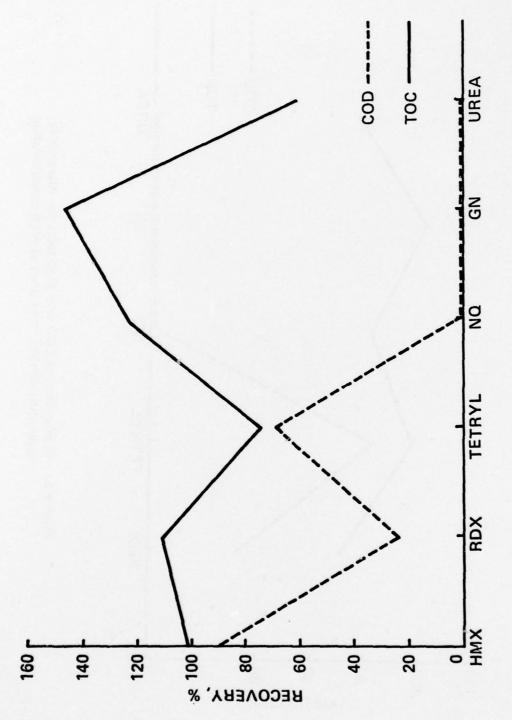


Figure 6. Comparison of COD and TOC tests for nitramines, guanidine nitrate, and urea at low concentration

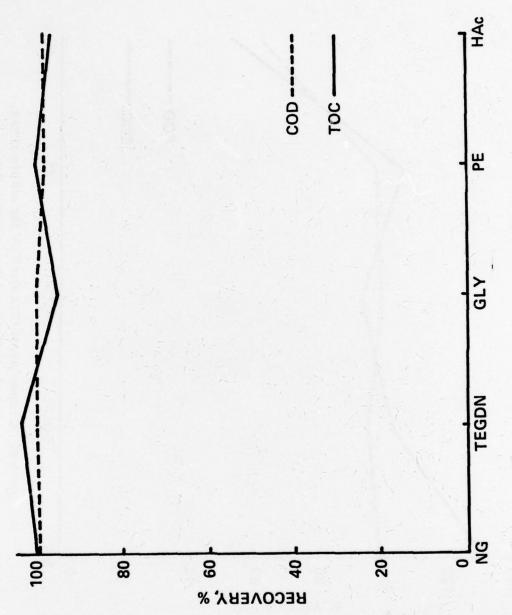

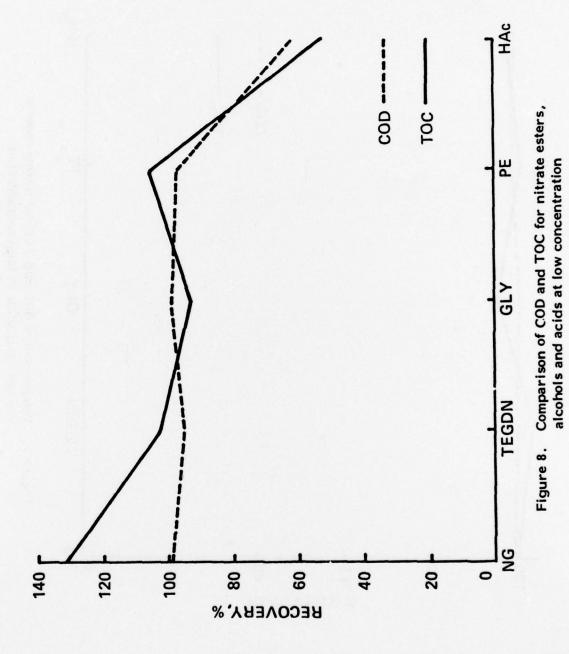

# **DISTRIBUTION LIST**

## Commander

US Army Armament Research and Development Command

ATTN: DRDAR-CG

DRDAR-LCM

DRDAR-LCM-SA (6)

DRDAR-LCP-F

DRDAR-TSS (5)

Dover, NJ 07801

# Commander

US Army Materiel Development

and Readiness Command

ATTN: DRCDE

DRCDE-E

DRCIS

DRCPA-E

DRCRP-I

DRCDE

DRCDL

DRCSG-S

5001 Eisenhower Avenue

Alexandria, VA 22333

# Commander

USDRC Installation & Services Agency

ATTN: DRCIS-RI-IU

DRCIS-RI-IC

Rock Island, IL 61299

### Commander

**US Army Armament Materiel** 

Readiness Command

ATTN: DRSAR-IMB-C (2)

DRSAR-RD

DRSAR-ISE (2)

DRSAR-SC

DRSAR-EN

DRSAR-PPW

DRSAR-ASF (2)

Rock Island, IL 61201

Project Manager for Munition Production
Base Modernization and Expansion
US Army Materiel Development
and Readiness Command
ATTN: DRCPM-PBM-EC
DRCPM-PBM-T-EV
Dover, NJ 07801

Department of the Army
Office, Chief of Research, Development
and Acquisition
ATTN: DAMA-CSM-P
Washington, DC 20310

Commander
US Army Procurement Equipment Agency
ATTN: DRX-PE-MT
Rock Island, IL 61299

Department of the Army ATTN: DAEN-ZCE Washington, DC 20310

Commander
Chemical System Laboratory
ARRADCOM
ATTN: DRDAR-CLM-T
Aberdeen Proving Ground, MD 21010

Defense Contract Administration Services (2) 1610 S. Federal Building 100 Liberty Avenue Pittsburgh, PA 15222

Defense Documentation Center (12) Cameron Station Alexandria, VA 22314 Commander
US Army Construction Engineering
Research Laboratory
ATTN: CERL-ER
Champaign, IL 61820

Office, Chief of Engineers ATTN: DAEN-MCZ-E Washington, DC 20314

US Army Engineer District, New York ATTN: Construction District 26 Federal Plaza New York, NY 10007

US Army Engineer District, Baltimore ATTN: Construction Division P.O. Box 1715 Baltimore, MD 21203

US Army Engineer District, Norfolk ATTN: Construction Division 803 Front Street Norfolk, VA 23510

US Army Engineer District, Mobile ATTN: Construction Division P.O. Box 2288 Mobile, AL 36628

US Army Engineer District, Fort Worth ATTN: Construction Division P.O. Box 17300 Fort Worth, TX 76102

US Army Engineer District, Omaha ATTN: Construction Division 6014 USPO and Courthouse 215 North 17th Street Omaha, NE 68102

US Army Engineer District, Kansas City ATTN: Construction Division (2) 700 Federal Building Kansas City, MO 64106 US Army Engineer District, Sacramento ATTN: Construction Division 650 Capitol Mall Sacramento, CA 95814

US Army Engineer District, Huntsville ATTN: Construction Division P.O. Box 1600 West Station Huntsville, AL 35807

Commander
US Army Environmental Hygiene Agency
ATTN: HSE-E (2)
Aberdeen Proving Ground, MD 21010

Commander
Badger Army Ammunition Plant
ATTN: SARBA-CE
Baraboo, WI 53913

Commander
Cornhusker Army Ammunition Plant
ATTN: SARCO-E
Grand Island, NB 68801

Commander
Holston Army Ammunition Plant
ATTN: SARHO-E
Kingsport, TN 37662

Commander
Indiana Army Ammunition Plant
ATTN: SARIN-OR
Charlestown, IN 47111

Commander
Naval Weapons Support Center
ATTN: Code 5042, Mr. C. W. Gilliam
Crane, IN 47522

Commander
Iowa Army Ammunition Plant
ATTN: SARIO-A
Middletown, IA 52638

Commander
Joliet Army Ammunition Plant
ATTN: SARJO-SS-E
Joliet, IL 60436

Commander
Kansas Army Ammunition Plant
ATTN: SARKA-CE
Parsons, KS 67537

Commander
Lone Star Army Ammunition Plant
ATTN: SARLS-IE
Texarkana, TX 57701

Commander
Longhorn Army Ammunition Plant
ATTN: SARLO-O
Marshall, TX 75670

Commander
Louisiana Army Ammunition Plant
ATTN: SARLA-S
Shreveport, LA 71102

Commander
Milan Army Ammunition Plant
ATTN: SARMI-S
Milan, TN 38358

Commander
Newport Army Ammunition Plant
ATTN: SARNE-S
Newport, IN 47966

Commander
Pine Bluff Arsenal
ATTN: SARPB-ETA
Pine Bluff, AR 71601

Commander
Radford Army Ammunition Plant
ATTN: SARRA-IE
Radford, VA 24141

Commander
Ravenna Army Ammunition Plant
Ravenna, OH 44266

Commander
Sunflower Army Ammunition Plant
ATTN: SARSU-O
Lawrence, KS 66044

Commander
Volunteer Army Ammunition Plant
ATTN: SARVO-T
Chattanooga, TN 34701

Army Logistics Management Center Environmental Management ATTN: Mr. Otto Nauman (2) Fort Lee, VA 23801

Project Manager for Chemical Demilitarization and Installation Restoration ATTN: DRCPM-DRR, Mr. Harry Sholk Aberdeen Proving Ground, MD 21010

Department of the Army ATTN: DAEN-FEU Washington, DC 20314

Dr. John A. Brown, Chairman P.O. Box 145 Berkeley Heights, NJ 07922

Dr. Helmut Wolf, Consultant 120 Skyline Drive Fayettesville, AR 72701

Dr. Fred Smetana, Consultant 5452 Parkwood Drive Raleight, NC 27612

Dr. Zachary Shetman, Consultant 109 N. Broadway White Plains, NY 10603 Commander
US Army Armament Materiel and
Readiness Command
ATTN: DRSAR-LEP-L
Rock Island, IL 61299

Director
US Army TRADOC Systems
Analysis Activity
ATTN: ATAA-SL (Tech Lib)
White Sands Missile Range, NM 88002

Weapon System Concept Team/CSL ATTN: DRDAR-ACW Aberdeen Proving Ground, MD 21010

Technical Library
ATTN: DRDAR-CLJ-L
Aberdeen Proving Ground, MD 21010

Technical Library
ATTN: DRDAR-TSB-S
Aberdeen Proving Ground, MD 21005

Technical Library
ATTN: DRDAR-LCB-TL
Benet Weapons Laboratory
Watervliet, NY 12189

# ED Records to the second of t